


Specimen Collection in Powerchart


1. Click on Power Chart on App Bar.




2. Click on Collection Runs.

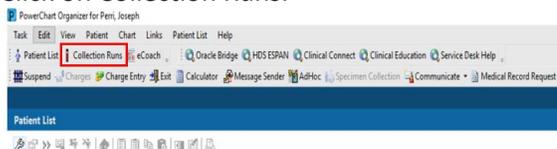


3. Add your locations using the search bar and ensure lab collect radial button is marked off.
4. Select specific clinical unit
5. Click View Collections icon to view patients.

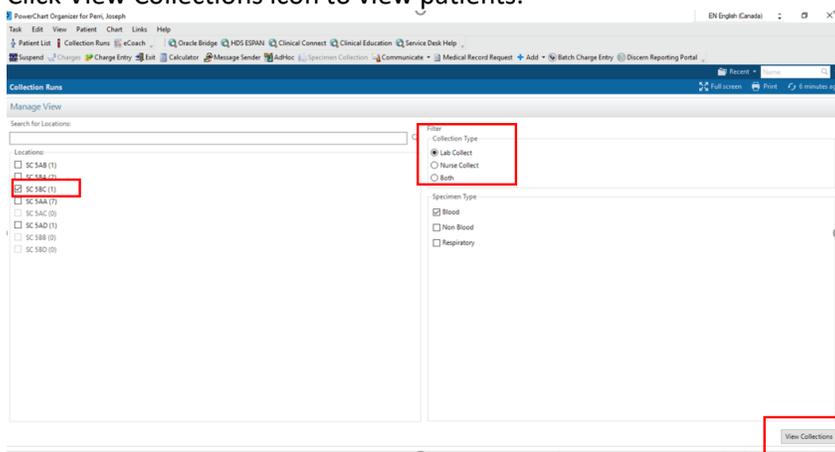

6. Click on Chart Collections.

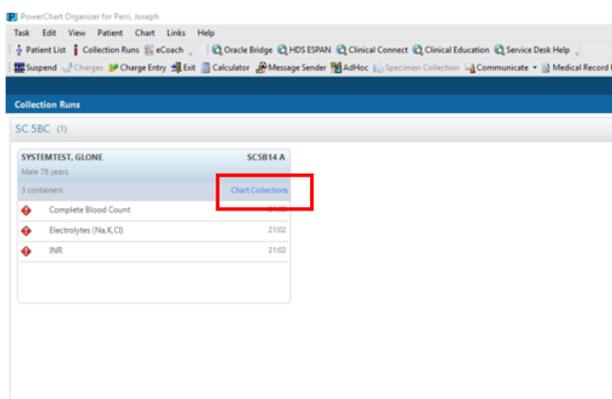


7. Scan Patient Barcode.

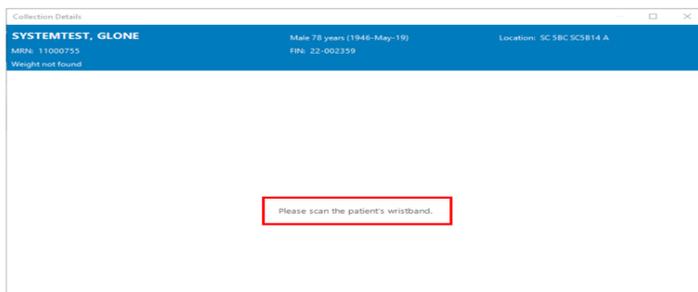


8. Perform patient ID and review specimens to be collected.

9. Print **all** collection labels using the printer icon or print **individual** labels using the pull down option for each specimen to be collected.

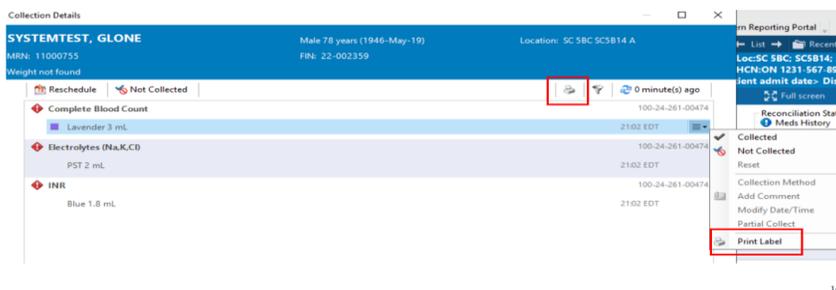

10. Mark each specimen as collected by scanning the specimen barcode or by using the pull down option for each specimen.

11. Click the sign icon to electronically sign and time stamp your collections.

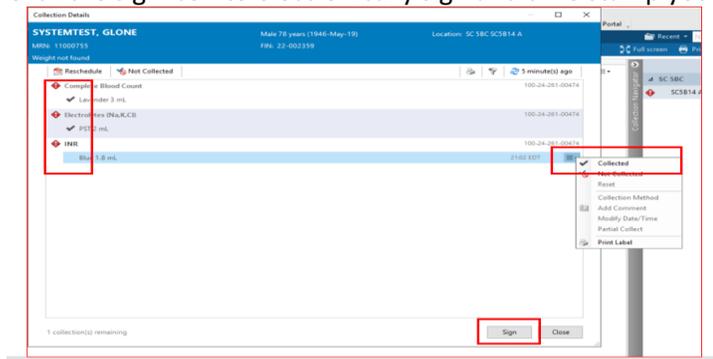

12. Transport specimens to lab for testing.